anti-GM1 autoantibodies

with Enzyme Labels IgG and IgM

ELISA

EK-GM1-GM 96 wells
(IgG & IgM)

Revision date: 2010-08-11

ENGLISH

INTENDED USE

BÜHLMANN anti-GM1 autoantibodies ELISA is designed for the quantitative determination of IgG and / or IgM isotypes of autoantibodies directed against GM1 in human serum (1-8) using individual IgG- and IgM-conjugates (2 results / patient).

PRINCIPLE OF THE ASSAY

BÜHLMANN anti-GM1 autoantibodies ELISA is based on the enzyme-immunometric assay technique. GM1 has been coated onto the wells of the microtiter plate (

Figure 1). Calibrator, controls and patient sera are incubated for two hours in the microtiter wells and anti-GM1 autoantibodies (Ab) present in the sample bind to the immobilized GM1. After washing off unbound substances, horseradish-peroxidase (HRP) labeled antibodies against human IgG and / or IgM are added to the wells and incubated for another two hours. Following a second washing step in which unbound antibody-enzyme reagent is removed, a substrate solution containing tetramethylbenzidine (TMB) is added to the wells. A blue color develops in proportion to the amount of anti-GM1 autoantibodies bound to the microplate in the initial step. Color development is stopped by adding an acidic stop solution which turns the blue solution into yellow. The intensity of the color is measured at 450 nm.

The measured absorbance is proportional to the titer of anti-ganglioside antibodies present in a given sample. The titers of anti-ganglioside autoantibodies are expressed as % ratios of the calibrator and can be assigned to titer categories.

REAGENTS SUPPLIED AND PREPARATION

Reagents	Quantity	Code	Reconstitution
Microtiter Plate precoated with GM1	12 x 8 wells	B-GM1-MP	Ready to use
Plate Sealer	3 pieces		
Wash Buffer Concentrate (10X) with preservatives	1 bottle 100 ml	B-GCO-WB	Dilute with 900 ml of deionized water
Incubation Buffer with preservatives	1 bottle 100 ml	B-GCO-IB	Ready to use
Calibrator Lyophilized with preservatives	1 vial	B-GCO-CA	Add 1.5 ml of Incubation Buffer
Negative, Low and Medium Control Lyophilized with preservatives	3 vials	B-GCO-CONSET	Add 1.5 ml of Incubation Buffer
Enzyme Label IgG Anti-human IgG Ab conjugated to HRP in a protein-base buffer with preservatives	1 vial 11 ml	B-GCO-ELG	Ready to use
Enzyme Label IgM Anti-human IgM Ab conjugated to HRP in a protein-base buffer with preservatives	1 vial 11 ml	B-GCO-ELM	Ready to use
TMB Substrate TMB in citrate buffer with H ₂ O ₂	1 vial 11 ml	B-TMB	Ready to use
Stop Solution 0.25 M sulfuric acid	1 vial 11 ml	B-ST5	Ready to use Corrosive agent

Table 1

STORAGE AND SHELF LIFE OF REAGENTS

Sealed / Unopened Reagents	
All sealed/unopened kit components are stable at 2-8°C until the expiration date printed on the labels.	
Opened / Reconstituted Reagents	
Microtiter Plate	Return unused strips immediately to the aluminium pouch containing the desiccant packs and reseal along the entire edge of zip-seal. Store for up to 4 months at 2-8°C.
Diluted Wash Buffer	Store for up to 4 months at 2-8°C.
Calibrator	Store for up to 1 month at 2-8°C. Do not freeze!
Controls	
Incubation Buffer	Store at 2-8°C until expiration date printed on the labels.
Enzyme Label	
TMB Substrate	
Stop Solution	Store at 18-28°C.

Table 2

WARNINGS AND PRECAUTIONS

Both, Calibrator (B-GCO-CA) and Controls (B-GCO-CONSET) of this kit contain components of human origin. Each serum donor unit used in the preparation of the components was tested by an FDA approved method and found negative for HBV surface antigen, so as for HCV and HIV1/2 antibodies. Although these methods are highly accurate, there is no guarantee that this material cannot transmit Hepatitis or AIDS. *Therefore, all patient specimens and kit components should be handled as if capable of transmitting infections.* All products containing human source material should be handled in accordance with good laboratory practice and using appropriate precautions.

MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipettes with disposable tips: 20 µl, 100 µl and 1 ml pipettes.
- Disposable polystyrene or polypropylene tubes for the preparation of sample dilutions.
- 1000 ml cylinder for the reconstitution of the wash buffer.
- Squeeze bottle for wash buffer or automatic microtiter plate washer.
- Blotting paper.
- Microtiter plate rotator.
- Microtiter plate reader for measurement of absorbance at 450 nm.

SPECIMEN COLLECTION AND STORAGE

- The procedure requires <0.1 ml of blood and <50 µl of serum, respectively.
- Refer to page 4 to learn about the interference of hemolyzed, lipemic or icteric samples.
- Collect blood into plain tubes (no anti-coagulant), avoid hemolysis, leave to clot for one hour, centrifuge for 10 minutes at approximately 1500 x g at room temperature (18-28°C), collect the serum.
- We recommend freezing aliquots of patient samples if you need to store samples in order to avoid repeated freezing/thawing.
- Store serum samples at ≤ -20°C up to 4 months. For long-term storage we recommend -70°C (samples are stable for >1 year). Frozen samples should be thawed and vortexed thoroughly prior to use.

PROCEDURAL NOTES

Recommendations and caveat for performance

Results and outcome of the assay may be influenced by many factors. It is important to strictly follow the protocol and the pre-analytical requirements. In particular we recommend attending to the following steps:

- **Residues in the microtiter plate wells** result from the production process. They are removed in the washing step (Assay procedure step 3) and do not affect the results.
- **Steps 3-9:** Use cold (2-8°C) reagents for all these steps and keep them cold while pipetting.
- **Steps 3, 6, 9:** Make sure that the wells are completely empty after the last washing cycle.
- **Step 9:** Adjust TMB Substrate to room temperature (18-28°C) before using it.
- **Step 11:** Shake microtiter plates during the incubation with substrate. Depending on the plate shaker, we recommend 400-600 rpm. The solution should be moved in the wells but must not spill over.
- If an **automated washer is used**, "plate mode" should be chosen so that dispensing is performed sequentially on all strips before aspirating.

ASSAY PROCEDURE

1. **Dilute** all patient samples to be investigated 1:50 with Incubation Buffer. Use 30 µl of serum + 1470 µl of Incubation Buffer. Mix by vortexing and leave diluted samples and reconstituted calibrator and controls for 30 minutes at 2-8°C prior to pipetting (refer to step 4 a-c).
2. **Prepare a plate-frame** with the required number of strips to test the patient samples. Reseal the remaining strips in the foil pouch together with the desiccant packs **immediately after usage**. Store refrigerated.

Note: Use cold reagents in steps 3 to 9.

3. **Wash** the coated wells twice using at least 300 µl of **cold** Wash Buffer per well. Empty the wells and tap the plate firmly onto blotting paper to remove remaining liquid completely.

Note: Immediately proceed to the next step (4a-e).

Detection of IgG-Isotype:

Note: We recommend testing Calibrator, Controls and samples in duplicates

- 4a. Calibrator: Pipet 100 µl of the **Calibrator** into the well A1 and A2 (see Figure 1).
- 4b. Controls: Pipet 100 µl of the **Medium Control** into well B1, and B2, **Low Control** into well C1 and C2 and **Negative Control** into the well D1 and D2 (see Figure 1).
- 4c. Patient serum: Pipet 100 µl of **diluted patient serum 1** into the wells E1 – E2 (see Figure 1).
- 4d. Patient serum: Pipet 100 µl of **diluted patient serum 2** into the wells F1 – F2 (see Figure 1).
- 4e. Pipet 100 µl of **diluted patient sera x-y** into the subsequent wells (see Figure 1).

Detection of IgM-Isotype:

5. Repeat step 4a – 4e using the subsequent wells.
6. Cover the plate with a Plate Sealer and **incubate** for 2 hours ± 5 minutes at 2-8°C (do not shake the plate).

7. Remove the Plate Sealer. Empty the wells and **wash** three times using at least 300 µl of **cold** Wash Buffer (2-8°C) per well. Empty the wells and strike the plate firmly onto blotting paper in order to remove washing buffer completely.
8. Add 100 µl of **Enzyme Label IgG or IgM** to the respective wells.
9. Cover the plate with a Plate Sealer and **incubate** for 2 hours ± 5 minutes at 2-8°C (do not shake the plate).
10. Remove the Plate Sealer. Empty the wells and **wash** three times using at least 300 µl of **cold** Wash Buffer (2-8°C) per well. Empty the wells and strike the plate firmly onto blotting paper.

Note: Let TMB Substrate Solution adjust to room temperature (18-28°C).

11. Add 100 µl of **TMB Substrate Solution** to each well.
12. Cover the plate with a Plate Sealer, incubate the plate on a **plate shaker** at 400-600 rpm for 30 minutes ± 2 minutes at 18-28°C protecting the plate from direct light.
13. Add 100 µl of **Stop Solution** to all wells. Proceed to step 14. within 30 minutes.
14. Read the **absorbance** at 450 nm in a microtiter plate reader.

RESULTS AND CALCULATION

Standardization

The Calibrator included in this kit has been calibrated against internal reference material. It has been **adjusted to 100 % ratio**.

Calculation of results:

1. Record the absorbance (OD) at 450 nm for each well (Calibrator, Controls and patient samples).
2. Average the duplicate Calibrator and Control values (if available).
3. Results are expressed as ratio of the absorbance of samples and the (averaged) absorbance of the Calibrator

$$\% \text{ Ratio} = \frac{\text{absorbance of samples and Controls}}{\text{absorbance of Calibrator}} \times 100$$

Programs to calculate results as % ratio are available on most microplate readers.

Note: The results presented in Table 11 are examples. Calibrator and Controls must be used in each individual assay.

QUALITY CONTROL

A good understanding of this instruction for use is necessary to obtain reliable results. These will be obtained only by using precise laboratory techniques (current GLP guidelines) and accurately following the instruction for use.

Since there is no control serum for anti-ganglioside antibodies commercially available, we recommend using a positive, and negative serum pool for internal quality control.

All controls must be within established confidence ranges (% ratio). The confidence ranges of the Controls are lot-specific and printed on the QC data sheet delivered with this kit.

Performance characteristics should be within established limits. If these characteristics are not in conformity with established limits and repetition excludes handling failures, check the following issues: i) Have all reagents, used in step 3-9, been kept at 2-8°C? ii) accuracy of the pipets, thermometers, and timers, iii) settings of ELISA washer and reader, iv) expiration date of the reagents v) storage and incubation conditions vi) color of the TMB Substrate Solution (should be colorless) vii) purity of the water.

PERFORMANCE CHARACTERISTICS

Intra-Assay Precision (Within-Run): 7.0 %. The intra-assay precision was calculated from results of 12 values of three IgM and IgG samples in a single run. The results in Table 12 have been calculated as % ratio as described in "Results and Calculation".

Inter-Assay Precision (Run-to-Run): 11.1 %. The inter-assay precision has been determined by measuring three serum samples with IgG and three with IgM antibodies to GM1 in 20 different runs. The results in

Table 13 have been calculated as % ratio as described in "Results and Calculation".

Analytical Sensitivity: 12 Incubation Buffer replicates were assayed in a single run. The detection limit expressed as the % ratio of the calibrator was calculated to be ≤ 5 %.

Linearity: The linear range of the test system was assessed according to CLSI guideline EP06-A. The system is linear in the diagnostic relevant range between 20 and 100 % ratio. Results above 100 % ratio are assessed clinically correct and can be diluted into the linear range by an additional 1:5/1:10 dilution.

Specificity: Different human serum samples containing specific anti-ganglioside IgM and/or IgG autoantibodies were incubated over night with the corresponding soluble antigen in different concentrations and subsequently tested in the BÜHLMANN Anti-GM1 Autoantikörper ELISA according to the assay procedure. Specificity of the autoantibody binding was demonstrated by inhibition with the corresponding antigen at concentrations between 1 and 100 µg/ml (data not shown).

INTERFERING SUBSTANCES

No interference is detected with the following substances up to the following concentrations: Triglycerides (Intralipid®): 3000 mg/dL; conjugated bilirubin: 60mg/dL; unconjugated bilirubin: 40 mg/dL and hemoglobin: 400 mg/dL..

EXPECTED VALUES AND CUT-OFF

In co-operation with the institutions mentioned below, a cut-off of 50% has been established. Values < 30% have been clearly classified as negative. The established titer categories are based on n=100 blood donors¹ (adult men and women between 18 to 70 years of age) and n=277 pathological samples². Sera were assayed for anti-GM1 autoantibodies according to the assay procedure. Results of normal blood donors are shown in Table 14.

Guidelines for the Use of Cut-off and Titer categories/ratios (%):

Category	Ratio (%)
Negative	< 30 %
Grey zone	30-50 %
Cut-off	50 %
Positive	> 50-100 %
Strongly Positive	> 100 %

¹) collected at the blood donation centre, University Hospital of Basel

²) received from Friedrich-Baur-Institute, Ludwig-Maximilians-University of Munich; Department of Neurology, University of Basel; Department of Neurology, University of Lyon.

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Table 14 **Blood Donors**

Category	anti-GM1-autoantibodies (%)	
	IgM	IgG
Negative	95	96
Greyzone	3	2
Positive	2	2
Strongly positive	0	0

Legend table 14: Distribution frequency of anti-GM1 autoantibodies (%) tested with EK-GM1-GM and classified into titer categories (n=100); calculated for IgG and IgM isotype:

- Negative : 95-96 % of blood donors
- Grey zone: 2-3 % of blood donors
- Positive: 2 % of blood donors
- Strongly positive: 0 % of blood donors

Table 12 **Intra-Assay Precision (Within-Run)**

Anti-GM1	Serum	Mean [% Ratio]	SD [% Ratio]	CV [%]
Enzyme-Label IgG	1	116	11.1	9.6
	2	60	5.0	8.3
	3	150	11.2	7.5
Enzyme-Label IgM	1	78	5.3	6.8
	2	60	2.0	3.4
	3	105	6.6	6.3
Mean IgG/IgM				7.0

Table description: cf. "Results and Calculations", "Performance Characteristics", and "Expected Values and Cut-Off" (pages 3/4)

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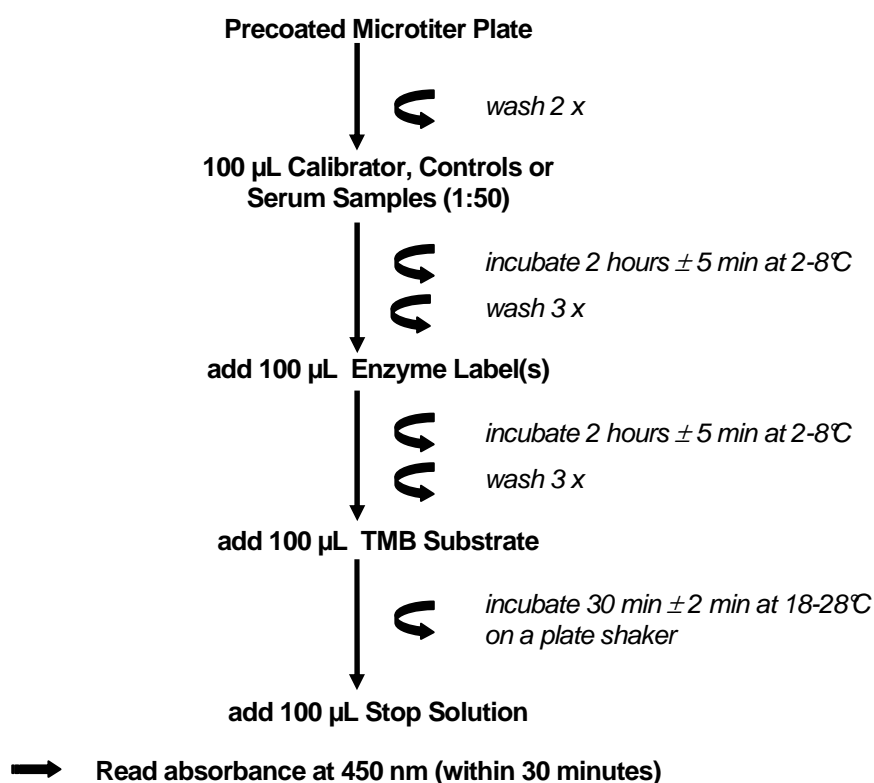
1. Willison HJ and Yuki N: *Peripheral neuropathies and anti-glycolipid antibodies*. Brain 125, 2591-2625 (2002).
2. Latov N: *Antibodies to glycoconjugates in neurological disease*. Clin Aspects Autoimm 4, 18-29 (1990).
3. Pestronk A: Invited Review: *Motor neuropathies, motor neuron disorders, and antiglycolipid antibodies*. Muscle Nerve 14, 927-936 (1991).
4. Steck AJ: *Use of anti-glycoconjugate antibody assays in neuropathy*. Workshop presentation at the Nat. Meeting Am. Acad. Neurology, Seattle, May 6-13, 1995.
5. Steck AJ and Kappos L: *Changing concepts in inflammatory and paraneoplastic neuropathies*. Curr Neurol 16, 191-212 (1996).
6. Humbel RL and Schmit P: *Anticorps antigangliosides et neuropathies peripheriques*. Rev Med Liege 51, 368-375 (1996).
7. Quarles RH and Weiss MD: *Autoantibodies associated with peripheral neuropathy*. Muscle Nerve 22, 800-822 (1999).
8. Nobile-Orazio E: et al: *How useful are anti-neural IgM antibodies in the diagnosis of chronic immune-mediated neuropathies?* JNS 266, 156-163 (2008).

BÜHLMANN anti-GM1 autoantibodies ELISA

TIME TO RESULT: 4.5 HOURS

SYMBOLS/ SYMBOLE/ SYMBOLES/ SIMBOLI/ SIMBOLOS

Symbol	Explanation
	Use By Verwendbar bis Utiliser jusqu'au Utilizzare entro Fecha de caducidad
REF	Catalogue number Bestellnummer Référence du catalogue Numero di catalogo Número de catálogo
LOT	Batch code Chargenbezeichnung Code du lot Codice del lotto Codigo de lote
IVD	<i>In Vitro</i> Diagnostic Medical Device <i>In Vitro</i> Diagnostikum Dispositif médical de diagnostic <i>in vitro</i> Dispositivo medico-diagnostico <i>in vitro</i> Producto sanitario para diagnóstico <i>in vitro</i>
	Temperature limitation Zulässiger Temperaturbereich Limites de température Limiti di temperatura Limite de temperatura
MP	Microtiterplate Mikrotiterplatte Plaque de microtitration Microplaca Microplaca
BUF WASH 10X	Wash Buffer Concentrate (10x) Waschpuffer-konzentrat (10x) Concentré de tampon de lavage Tampón de lavado concentrado (x10) Tampón de lavado concentrado (x10)
BUF INC	Incubation Buffer Inkubations-Puffer Tampon d'incubation Tampón de incubación Tampón de incubación

Symbol	Explanation
CONTROL-	Negative Control Negativkontrolle Contrôle négatif Control negativo Control negativo
CONTROL L	Low Control Kontrolle tief Contrôle faible Control bajo Control bajo
CONTROL M	Medium Control Kontrolle mittel Contrôle moyen Control medio Control medio
CAL	Calibrator Kalibrator Calibrateur Calibratore Calibrador
EL IgG	Enzyme Label IgG Enzymmarker-IgG Marqueur enzymatique IgG Marcato enzimatico IgG Marcador enzimático IgG
EL IgM	Enzyme Label IgM Enzymmarker-IgM Marqueur enzymatique IgM Marcato enzimatico IgM Marcador enzimático IgM
SUBS TMB	TMB Substrate TMB Substrat Substrat TMB Substrato di TMB Substrato de TMB
SOLN STOP	Stop Solution Stopp-Lösung Solution stop Soluzione stoppante Solución de parada

Printing Date
2010-11-17