

Quantum Blue® Adalimumab

trough level rapid test

Modulare
la terapia
senza ritardi



- Il primo test rapido quantitativo per il monitoraggio di Adalimumab
- Risultato in 15 minuti
- Range 1,3 – 35 µg/mL
- Elevata concordanza con i test ELISA di riferimento


Quantum Blue® Adalimumab test rapido

Caratteristiche principali

- Il primo test rapido quantitativo per Adalimumab
- Range di misura 1,3 – 35 µg/mL
- Copre il range di 5 – 12 µg/mL rilevante per il Trough Level*
- Ottima correlazione con i test ELISA di riferimento

* Mitrev et al. 2017

Correlazione con test ELISA:

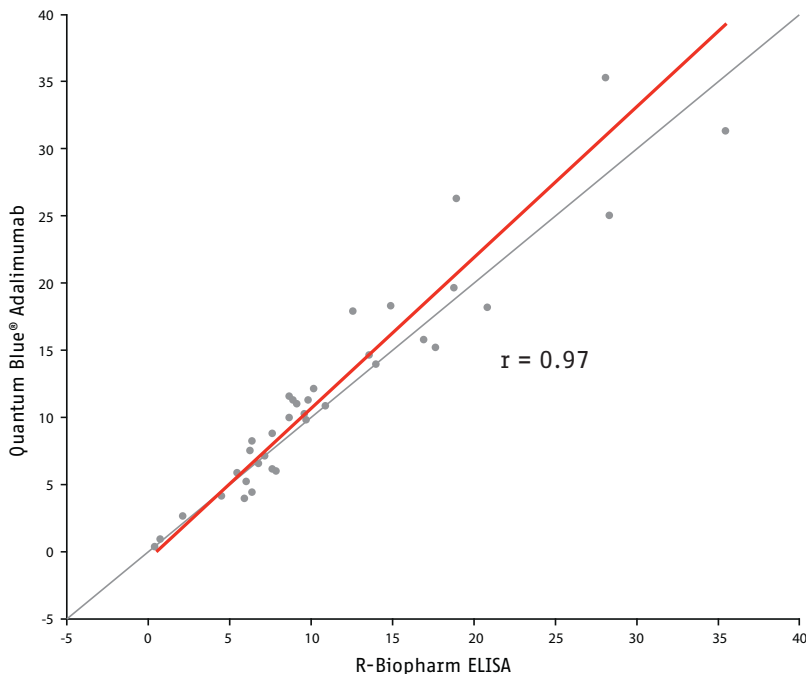

— Passing-Bablok fit

		RIDASCREEN ADM MONITORING (ELISA, µg/mL)		
		low	optimal	high
Quantum Blue® Adalimumab (µg/mL)	low	20.0%	5.0%	0.0%
	optimal	0.0%	27.5%	0.0%
	high	0.0%	12.5%	35.0%

Concordanza totale 82.5%



BÜHLMANN Laboratories AG
Switzerland
info@buhlmannlabs.ch
www.buhlmannlabs.ch
Phone +41 61 487 12 12

BÜHLMANN Italia Srl
Italy
info@buhlmannlabs.it
www.buhlmannlabs.it
Fax ordini +41 61 487 12 99

Quantum Blue® is a registered trademark
of BÜHLMANN in many countries.

Codici ordini:
LF-TLAD10 10 card
LF-TLAD25 25 card

Prodotti correlati:
LF-TLIF10 10 card
LF-TLIF25 25 card
LF-CAL25 25 card
LF-CALE25 25 card
LF-CHR25 25 card

CE-marked products